


# CELLCOAT®

## Poly-L-Lysine coated cell culture vessels

Poly-D-Lysine and Poly-L-Lysine are synthetic molecules which increase cell adhesion to vessels by altering the net charge of the surface without stimulating biological activity. Applications include cultivation of transfected cells and other fastidious cell lines, serum-free or –reduced cell culture and general improvement of cell attachment e.g. during washing steps.

Furthermore PLL improve survival of many neuronal cells in culture. As synthetic molecule PLL is free of impurities carried by other proteins.

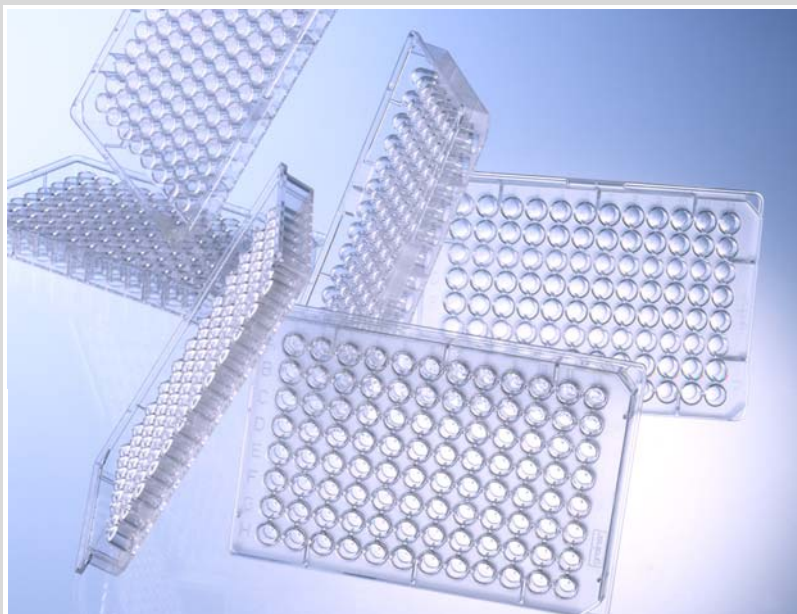

### Features

Source	synthetic
MW	PDL : 75-150 kD PLL : 30-70 kD
QA	<ul style="list-style-type: none"><li>■ Tested for ability to promote attachment of a neuroblastoma cell line</li><li>■ Negative for bacteria and fungi</li></ul>
Applications	<ul style="list-style-type: none"><li>■ Cultivation of transfected cells and other fastidious cell lines</li><li>■ Serum-free and serum-reduced culture</li><li>■ Improvement of cell attachment</li></ul>

### Ordering Information

Cat.-No.	Poly-L-Lysine coated cell culture vessels	Qty/bag	Qty/case
628 930	60 mm dish, PLL coated	20	100
657 930	6 well multiwell plate, PLL coated	5	50
662 930	24 well multiwell plate, PLL coated	5	50
655 930	96 well microplate, clear, PLL coated	5	20
655 936	96 well microplate, µclear®, black, PLL coated	5	20
781 930	384 well microplate, clear, µClear PLL coated	5	20
781 936	384 well microplate, black, µClear® PLL coated	5	20

Revision: October, 2005 – 075 016

[www.gbo.com/bioscience](http://www.gbo.com/bioscience)

**Germany (Main office)**  
Greiner Bio-One GmbH  
Phone: (+49) 7022 948-0  
E-Mail: [info@de.gbo.com](mailto:info@de.gbo.com)

**Belgium**  
Greiner Bio-One N. V.  
Phone: (+32) 2-4 61 09 10  
E-Mail: [info@be.gbo.com](mailto:info@be.gbo.com)

**France**  
Greiner Bio-One SAS  
Phone: (+33) 1 69-86 25 50  
E-Mail: [infos@fr.gbo.com](mailto:infos@fr.gbo.com)

**Japan**  
Greiner Bio-One Co. Ltd.  
Phone: (+81) 3-35 05-88 75  
E-Mail: [info@jp.gbo.com](mailto:info@jp.gbo.com)

**Austria**  
Greiner Bio-One GmbH  
Phone: (+43) 7583 6791-0  
E-Mail: [office@at.gbo.com](mailto:office@at.gbo.com)

**Netherlands**  
Greiner Bio-One B. V.  
Phone: (+31) 172-42 09 00  
E-Mail: [info@nl.gbo.com](mailto:info@nl.gbo.com)

**UK**  
Greiner Bio-One Ltd.  
Phone: (+44) 1453-82 52 55  
E-Mail: [info@uk.gbo.com](mailto:info@uk.gbo.com)

**USA**  
Greiner Bio-One Inc.  
Phone: (+1) 704-261-78 00  
E-Mail: [info@us.gbo.com](mailto:info@us.gbo.com)



**greiner bio-one**